



CliniWax Sutures

Description

CliniWax is a sterile mixture of beeswax, paraffin, and isopropyl palmitate.

Intended use

CliniWax sutures are intended for use in controlling the bleeding from bone surfaces.

Indications

CliniWax is indicated to stop bleeding from bones and it does not contain any medicine to stop bleeding

Actions

CliniWax acts as a mechanical barrier to achieve local hemistasis of the bone.

Contra-indications

The product should not be used where rapid osseous regeneration and fusion are desired. CliniWax should not be used on patients with known sensitivities or allergies to beeswax, paraffin wax and iso propyl palmitate.

Warning

1. Should not be exposed to excessive heat.
2. Do not reuse, reprocess or re-sterilize.
3. Should not be used if the packing is open or damaged.
4. Do not manipulate with latex gloves as latex proteins deposited in the Bone wax may precipitate an allergic reaction.

Precautions

CliniWax should be used sparingly. Excess CliniWax must be removed from the operating site. The package should be open prior to use to minimize the possibility of contamination and excessive drying. CliniWax melts at 62 °C – 64 °C.

Adverse Reactions

Adverse effects associated with the use of this device include mild inflammatory in tissues immediately adjacent to the site of implantation. CliniWax may enhance an existing infection and pollen granules imbedded in beeswax may elicit an allergic reaction.

Sterility

CliniWax sutures are sterilized by gamma radiation. Do not re-sterilize. Do not use if package is opened or damaged. Discard opened unused sutures.

Storage


Recommended storage conditions: between 15°C to 30°C, away from direct heat. Do not use after expiry date.

How supplied

CliniWax is supplied at a mass of 2.5 gram.



AKACIA MEDICAL (PTY) LTD
25 Harrower Road
North End
Port Elizabeth, 6001
Republic of South Africa
Tel: (+27) 41 373 9599
Email: clinisut@akaciamedical.co.za/
csdpe@akaciamedical.co.za
Website: www.clinisut.com
Made in South Africa


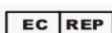


Obelis S.A.
Bd. Général Wahis 53
B-1030 Brussels, Belgium

Document No: PMF 12.10.3
Revision: 0
Date: 03 December 2024

CliniWax Sutures


INSTRUCTIONS FOR OPENING OF CLINISORB SUTURE

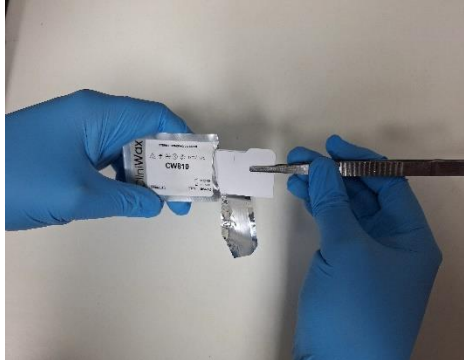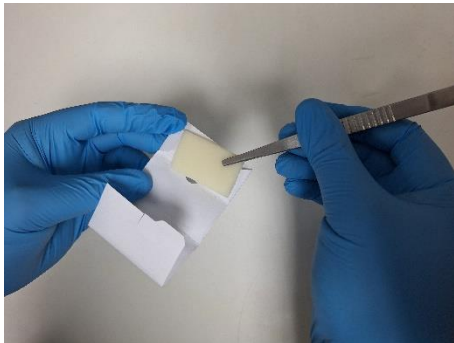1. Hold suture facing you with opening on right side. 	2. Open outer pouch by pulling it apart until inner pouch is exposed and remove. 
3. Open foil pouch at tearing notch. 	4. Tear open the foil pouch and pull the winding pouch out. 
5. Open the winding pouch. Secure the bone wax with tweezers from the winding pouch and pull it out. 	The suture is now ready for use






CliniWax Sutures

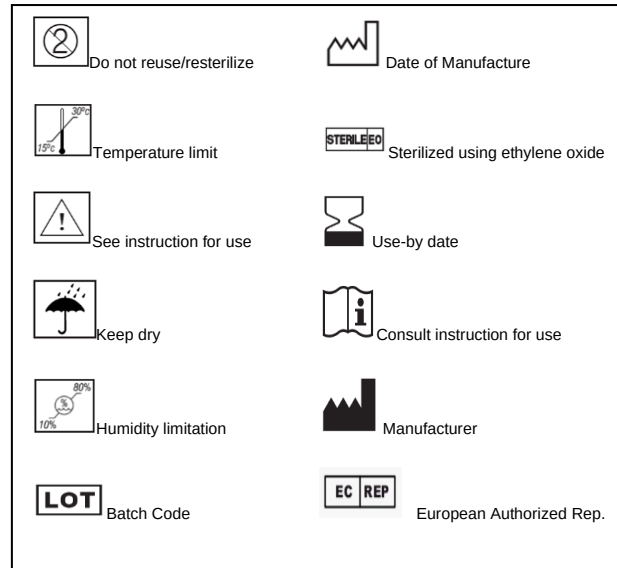

AKACIA MEDICAL (PTY) LTD
25 Harrower Road
North End
Port Elizabeth, 6001
Republic of South Africa
Tel: (+27) 41 373 9599
Email: clinisut@akaciamedical.co.za/
csdpe@akaciamedical.co.za
Website: www.clinisut.com
Made in South Africa

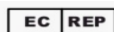

Obelis S.A.
Bd. Général Wahis 53
B-1030 Brussels, Belgium

Document No: PMF 12.10.3
Revision: 0
Date: 03 December 2024